



How the Baja Back Brace Line is Different

Baja – Posterior Compression

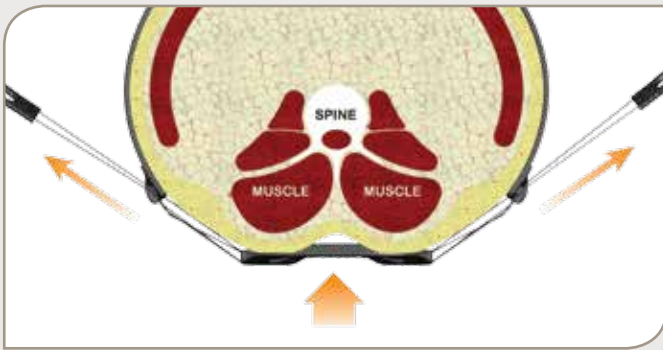

The Baja has a bi-lateral lacing system with two banks of pulleys positioned on the sides of the torso. When the patient engages the cinching mechanism, the posterior panel is drawn forward, focusing the compression on the soft tissues that support the spine.

Other Brands – Lateral Compression

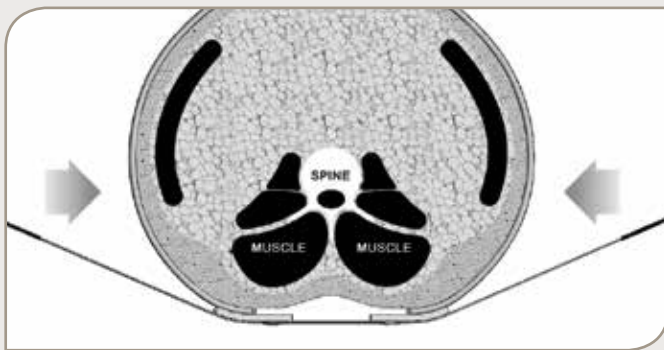

Other cinch-style back braces have a posterior lacing system with the laces positioned across the spine. When a patient engages this style of brace, the two sides pull together and the compression is focused on the sides of the torso.


Help your patients stay active.
Baja back braces provide immediate pain relief with targeted compression.



Relieve Back Pain
Restore Activity



Sizing Chart

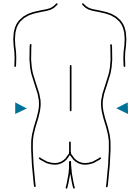	Waist measurement			
	Size	*PN	Inches	Centimeters
	S	1	26 – 30	66 – 76
	M	2	30 – 36	76 – 91
	L	3	36 – 44	91 – 112
	XL	4	44 – 54	112 – 137
	UNIV	4U	26 – 54	66 – 137

Products

Description	Product #
Baja 627 Back Brace <ul style="list-style-type: none">• 2 hot/cold gel pads• Breathable foam pad• Trim to fit - universal option	3750*
Baja 631 Back Brace <ul style="list-style-type: none">• 2 hot/cold gel pads• Breathable foam pad• Trim to fit - universal option	3770*
Baja 637 Back Brace <ul style="list-style-type: none">• 2 hot/cold gel pads• Breathable foam pad• Trim to fit - universal option	3760*

Accessories

Hot/cold gel pad	97322
Extension Panel	97363



240 East Hersey St., Suite 2
Ashland, OR 97520-5202
t: 800-541-2455
f: 541-482-2119
CropperMed.com



Patent #: 8,409,122
©2015 Cropper Medical, Inc.
07/2015 PN: 12340D



Be Better™

Baja™
627

Easy-cinching, targeted compression for back pain

The Baja 627 provides unrivaled lumbar support and pain relief. The bi-lateral pulley system is unique in its ability to provide vectored compression to the lower spine. As the pulleys are engaged, the lumbar panel is drawn in toward the spine, targeting the compression to the source of pain.

INDICATIONS

- Herniated Discs
- Degenerative Disc Disease
- Lumbar Sprains and Strains
- Spondylolisthesis
- Osteoarthritis
- Radicular Pain

HCPCS: L0627, L0642



Baja™
631

Immediate, conservative pain relief with added stabilization

The Baja 631 provides immediate pain relief and extra stability throughout the sacral, lumbar, and lower thoracic spine. The bi-lateral cinching mechanism and independent lumbar panel provide targeted compression for pain relief in the vertebral column and the surrounding soft tissues.

INDICATIONS

- Herniated Discs
- Degenerative Disc Disease
- Spondylolisthesis
- Osteoarthritis

HCPCS: L0631, L0648



Baja™
637

Pain relief and coronal control

The Baja 637 uses the same lumbar panel as the Baja 631 to provide pain relief and sagittal control but also includes rigid side panels to limit coronal motion. The side panels may be removed to step down patients' treatment over time.

INDICATIONS

- Herniated Discs
- Degenerative Disc Disease
- Spondylolisthesis
- Osteoarthritis
- Lumbar Spinal Stenosis
- Post-Surgical Stabilization

HCPCS: L0637, L0650



41% Reduction in Low-back Pain*

- Immediate back pain relief
- Targeted compression
- Hot & cold gel pads

Conservative Treatment

For patients with mild to moderate back pain, the Baja line provides immediate pain relief with targeted compression. It's the most conservative treatment to help them manage their short-term pain.

For patients dealing with chronic pain, the Baja line provides you with a non-pharmacological tool to manage their pain as you work together to find long-term solutions.

Clinically Proven Pain Relief

A recent study using the bi-lateral cinching mechanism of the Baja showed an average reduction in VAS scores of 40% and a 41% improvement in Roland Morris Disability Questionnaire scores in patients with chronic low back pain.*

Prospective Study of Lumbar Back Brace In an Interventional Pain Practice. Rosenthal, Richard M.D. & Spencer, Shawn, BS. Practical Pain Management, (July/August 2011)

Integrated Hot/Cold Therapy



Each Baja includes two hot/cold gel pads which can be heated in the microwave for soothing heat, or cooled in the freezer for pain-relieving cold. The gel pads slide directly into the lumbar pad.